




ISOLUX



IsoLED Mini with or without permanent attachment to glasses

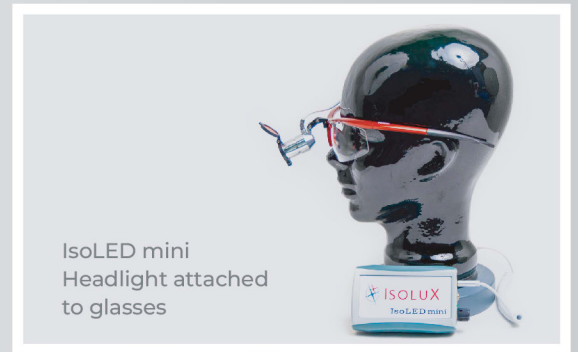

IsoLED mini Headlight attached to glasses



IsoLED mini Headlight including Safety Glasses

IsoLED Mini LED HEADLIGHT

The IsoLED Mini Portable LED Headlight is a Battery operated Exam / Diagnostic lamp ideal for use in the following environments:

- Dermatology
- Family Practice
- Veterinary
- Emergency Rooms
- Outpatient Facilities
- Urgent Care Centers
- Dentistry
- Minor Procedures

Its high efficiency design allows for reduced energy consumption and reduced cost of ownership (around 1/3 of halogen lamps), while its 40,000 hour of LED life provides for a free maintenance device with an end life of around 20 years.

The unit has infinite mobility, not confined by access to power sockets, and its charge will last more than 10 hours at full intensity.

Its light beam is extremely homogeneous and free of imperfections. The light- weight battery pack can be used on the user's belt or pocket. The Battery Charger comes with a Charge Monitor that alerts the operator when the Battery is fully charged.

The Headlight adapts to any Binocular loupes: Gallilei, Zumax, Zeiss, Orascoptic Flip up, Grendel Meier and Heine. The unit comes in a carrying case containing the Headlight with pendant cable/ connector, Battery Pack, two retaining multi-clips, removable yellow filter, Battery Charger, safety glasses and User's Guide.

239-514-7475 www.isoluxllc.com

1045 Collier Center Way Suite 6 & 7 Naples, FL 34110

SPECIFICATIONS

Illumination	> 22,000 lux @ 30 cm
LED Life	>40,000 hours
Color Temperature	4000 °K
Color Rendition Index (CRI)	>92Typical
Light Degradation	None
Spot Size	60 mm @ 300 mm working distance
Intensity Range	3% to 100%
Operating Time	>10 hours typical at maximum intensity
Normal Recharge Time	3.6 hours
Quick Recharge Time	2 hours, for 8 hours of maximum intensity
Battery Charger	Universal: 100 to 264 VAC, 50 to 60 Hz., 8.4 VDC @ 1.3 A
Weight (headlight)	23 grams
Weight (Battery Pack)	190 grams
Size (control unit/ battery pack)	4.25" x 2.7" x 1.1"
Warranty	2 Years on Headlight, 1 Year on Battery Pack